





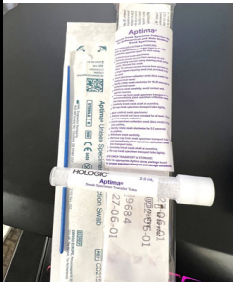








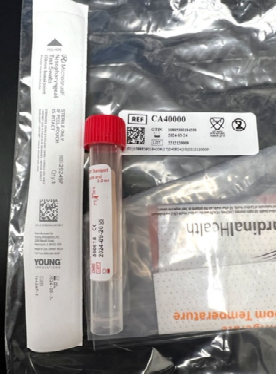
Device	Recommended Uses	Acceptable but Not Recommended	NEVER USED	Device name /Item Class & Stock #/Warehouse Description
Copan flocked E-swab 	Aerobic bacterial culture Anaerobic bacterial culture Yeast culture Gram Stain Fungal Stain Herpes & Varicella DNA probe (PCR) Viral PCR Respiratory syncytial Virus antigen Influenza A+B antigen testing	Skin Scrapings, send tissue Mycobacteria (AFB) cultures send unpreserved patient sample or tissue Fungal Culture send unpreserved patient sample or tissue Fluid Culture send unpreserved patient sample Tissue culture send unpreserved patient sample or tissue Operative samples send unpreserved patient sample or tissue Rapid Strep testing use EZ II culturette swab Respiratory Culture Wet Preparation	Bacteriology PCR Mycology PCR Viral Culture Stool Culture Chlamydia nucleic acid testing (DNA/RNA) Gonorrhea nucleic acid testing (DNA/RNA) HPV nucleic acid testing Trichomomas nucleic acid testing	Copan flocked E-swab Item Class: Misc. Patient Care Supplies: 554409 Description: SWAB COPAN ESWAB 480C
BD EZ II swab 	Rapid Strep testing Confirmatory streptococci culture COVID, Rapid NAAT Influenza A+B NAAT RSV NAAT	Virus culture, unless immediately placed in M5 buffer Viral PCR unless immediately placed in M5 buffer Mycoplasma PCR unless immediately placed in M5 buffer Ureoplasma PCR unless immediately placed in M5 buffer Aerobic bacterial Culture use Copan flocked E-swab Skin Scrapings, send tissue Mycobacteria (AFB) cultures send unpreserved patient sample or tissue Fungal Culture send unpreserved patient sample or tissue Fluid Culture send unpreserved patient sample Tissue culture send unpreserved patient sample or tissue Operative samples send unpreserved patient sample or tissue Respiratory syncytial Virus antigen Influenza A+B antigen testing	Anaerobic culture Bacteriology PCR Mycology PCR Stool Culture Chlamydia nucleic acid testing (DNA/RNA) Gonorrhea nucleic acid testing (DNA/RNA) HPV nucleic acid testing Trichomomas nucleic acid testing	BD EZ II swab Item Class: Misc. Patient Care Supplies: 554409 Description: SWAB BD EZII DUAL 220145
Sterile containers suitable for Aerobic / Anaerobic Culture 	Aerobic Culture Anaerobic culture Gram Stain Fungal Stain Mycobacteria (AFB) stain Mycobacteria (AFB) culture Fungal/Yeast Culture Fluid culture	Virus Culture send unpreserved patient sample	Viral PCR Bacteriology PCR Mycology PCR Wet preparation Cytology Sterile syringe <i>with</i> Needle	A.C.T II (Anaerobic transport) Item Class: Laboratory Supplies: 21608 Description: ACT II Dual Collection System 21608

Device	Recommended Uses	Not Recommended	NEVER USED	Device name /Item Class & Stock #/Warehouse Description
<p>Sterile container suitable for Aerobic Culture</p> 	<p>Aerobic Culture Urine Culture Sputum culture Gram Stain Fungal Stain Mycobacteria (AFB) culture Mycobacteria (AFB) stain Fungal/Yeast Culture Wet preparation Legionella antigen testing (Urine) Streptococcus pneumoniae antigen testing (Urine) Legionella culture Cytology Mycobacterial CX Legionella PCR</p>		<p>Anaerobic Culture Anaerobic Culture, ONLY IF tissue present</p>	<p>Item Class: Laboratory Supplies: 466001 Container Specimen 3 oz 466001 Item Class: Laboratory Supplies: 470106 LAB Centrifuge tube 25/pk 470106</p>
<p>Clean container</p> 	<p>Clostridium difficile Wet Preparation Ova and Parasite Other Fungal Culture</p>	<p>Aerobic culture, send unpreserved patient sample</p>	<p>Aerobic culture Anaerobic Culture Stool culture</p>	<p>Item Class: Laboratory Supplies: 466001 Container Specimen 3 oz 466001 Item Class: Laboratory Supplies: 470138/470139 Tube Lab Specimen Urine - 12 ml Cap Lab Spec tube 500/ca</p>
<p><i>Hologic Aptima</i> Urine Specimen Collection Kit Initial sample, not mid-stream</p> 	<p>Chlamydia trachomatis antigen Neisseria gonorrhoeae antigen "Nucleic Acid Amplification - RNA" In-house Chlamydia, GC testing for urine Mayo Trichomonas testing for urine</p>		<p>Sexual Abuse Rape Chlamydia Culture / PCR GC Culture / PCR Insufficient Sample Excess Sample HPV nucleic acid testing Trichomonas nucleic acid testing --> Mayo testing Aerobic / Anaerobic / Fungal / AFB culture</p>	<p><i>Hologic Aptima:</i> Urine Specimen Collection Kit Item Class: Laboratory Supplies: 51290 Description: GENPROBE URINE COLLECTION 1040</p>

Device	Recommended Uses	Not Recommended	NEVER USED	Device name /Item Class & Stock #/Warehouse Description
<p><i>Hologic Aptima</i> Unisex Specimen Collection Kit</p> <p>Endocervical and Male urethral samples</p> 	<p>Chlamydia trachomatis antigen Neisseria gonorrhoeae antigen</p> <p>"Nucleic Acid Amplification - RNA"</p> <p>Trichomonas nucleic acid testing</p> <p>In-house Chlamydia, GC and Trich testing for cervix and urethra</p>		<p>Sexual Abuse Rape Chlamydia Culture / PCR</p> <p>GC Culture / PCR White Swab</p> <p>White and Blue Swab</p> <p>HPV nucleic acid testing</p> <p>Aerobic / Anaerobic / Fungal / AFB culture</p>	<p><i>Hologic Aptima</i> Unisex Specimen Collection Kit</p> <p>Item Class: Laboratory Supplies: 51291</p> <p>Description: GENPROBE SPECIMEN SWAB 1041</p>
Device	Recommended Uses	Not Recommended	NEVER USED	Device name /Item Class & Stock #/Warehouse Description
<p>Hologic Aptima</p> <p>Multitest Specimen Collection Kit</p> 	<p>Chlamydia trachomatis antigen Oral or rectal Neisseria gonorrhoeae antigen Oral or Rectal ARUP Vaginitis Panel In-house Chlamydia and GC testing for Vaginal, Throat, and Rectal Swab Specimens</p>			<p>Item Class: Laboratory Supplies: 60627</p>

<p>Hologic - The "ThinPrep" Pap Test</p> 	<p>Chlamydia trachomatis antigen Neisseria gonorrhoeae antigen</p> <p>HPV antigen HPV genotyping HPV-16, 18/45</p> <p>Trichomonas nucleic acid testing</p>		<p>Sexual Abuse Rape Chlamydia Culture / PCR GC Culture / PCR</p> <p>Aerobic / Anaerobic / Fungal / AFB culture</p>	<p>Item Class: Laboratory Supplies: 70136 Thin Prep PAP test / 25 per pack</p>
<p>Alpha Tek ETM</p> 	<p>GI Diarrheal Panel</p> <p>Clostridium difficile</p> <p>Ova & Parasite: Cryptosporidium, Giardia, Cyclospora, Entamoeba histolytica</p>		<p>Ova and Parasite - Microscopic</p>	<p>Alpha Tek ETM</p> <p>Item Class: Laboratory Supplies: 900612</p> <p>Description: PARASET ETM ALPHATEC</p>

Device	Recommended Uses	Not Recommended	NEVER USED	Device name /Item Class & Stock #/Warehouse Description
	Ova and Parasite - Microscopic (Contact Microbiology at 3707 to place the order)		Clostridium difficile GI Diarrheal Panel Aerobic / Anaerobic / Fungal / AFB culture	Alpha Tek PROTO CLR Item Class: Laboratory Supplies: 400405 Description: PROTOFIX CLR 4602CS
Device	Recommended Uses	Not Recommended	NEVER USED	Device name /Item Class & Stock #/Warehouse Description
	Aerobic Anaerobic Blood Culture	Fluid Culture	Fungal Culture Mycobacteria culture	Item Class: Laboratory Supplies: 59268 BACTALERT Aerobic SA Blue Item Class: Laboratory Supplies: 59269 BACTALERT Anaerobic Pink

Device	Recommended Uses	Not Recommended	NEVER USED	Device name /Item Class & Stock #/Warehouse Description
	Fungal Blood Culture Mycobacteria Blood Culture		Aerobic Culture Anaerobic culture	Item Class: Laboratory Supplies: 51582 Vacutainer Green 6ml
Device	Recommended Uses	Not Recommended	NEVER USED	Device name /Item Class & Stock #/Warehouse Description
	COVID, Panther or Mayo COVID Flu Combo testing Respiratory Panel PCR, NP			Item Class: Laboratory Supplies